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Index

	page
Introduccion	21
Use	21
Where to find it	22
Caring for your Solar Charger	22
▲ Do not insert batteries with the incorret polarity	22
▲ Do not insert hearing aid with the incorret polarity	22
▲ Do not attempt to recharge batteries other than specified rechargeable batteries	23
▲ Only the following rechargeable batteries may be used with this charger	24
▲ Only use batteries with the codes described bellow or rechargeable batteries of equivalent size the hearing aid	24
▲ Do not disassemble	24
▲ Observe cleaning instructions	24
▲ Do not immerse in water	24
▲ Do not incinerate	25
▲ Place in direct sunlight to recharge	25
▲ If charging using household light	25
▲ Only government approved	25
▲ During charging the LED will illuminate	26
▲ Keep out of reach of small children, if batteries are swallowed seek medical attention immediately	26
Using your Solar Charger	26
Inserting Batteries	26
Inserting Solar Ear Rechargeable Hearing Aid	27
Removing Batteries	27
Warranty	28 e 29

Introduccion

The Solar Ear battery charger is designed to charge Solar Ear rechargeable hearing aids and rechargeable Ni-Mh hearing aids batteries #13, #312 and #675 as well as Ni-Mh AA batteries.

Standard zinc air hearing aid batteries can not be recharged. The Solar Charger can be used to recharge a hearing aid or two rechargeable batteries. Back-up batteries must be charged. If solar charger is opened and terminal components tampered with the warranty will be null and void.

It is recommended that the internal AA rechargeable batteries are charged in the day. Hearing aids and rechargeable batteries should be recharged at night.

Use

2 Internal AA Batteries

Not visible
(Charge for 20 hours
before first use)


Hearing Aid

Rechargeable
Hearing Aid
Batteries
(Flat side shows +)

Where to find it


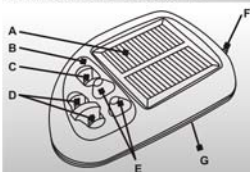
- A - Solar Panel
- B - Signalling LED
- C - Hearing Aid Recess
- D - Battery Recesses
- E - Battery Eject Buttons
- F - Input Electrical Charger
- G - 2 Internal AA Ni-Mh batteries

Caring for your Solar Charger

▲ Do not insert batteries with the incorret polarity.

Check the (+,-) on the charger when inserting the batteries. If the battery is inserted the wrong way around it may damage the charger. The LED will light up red if it happens.

The flat side of the battery is the + (positive) side.

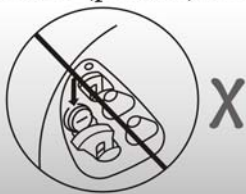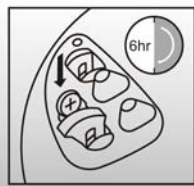
▲ Do not insert hearing aid with the incorret polarity.

If the hearing aid is inserted the wrong way around it may damage charger with reversed polarity can cause damage to the charger. The LED will light up red if it happen. Only Solar Ear hearing aids and Ni-Mh hearing aid or AA rechargeable batteries to be used with this charger.

Analogue

Digital

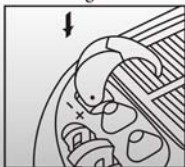If no contact is made, take the battery out of the hearing aid and charge it in the battery holder. For all other hearing aids which use #13, #312 or #675 rechargeable Ni-Mh batteries please use the battery holder for charging.

▲ Do not attempt to recharge batteries other than specified.

The manufacturer will not be held responsible should the incorrect batteries be charged with this charger.

▲ The following rechargeable batteries can be used with this charger.

Only a rechargeable #13, #312, #675 or equivalent rechargeable hearing aid battery can be used. Zinc Air, mercury, Ni-Cad or other button cell batteries will cause damage to the charger and must not be used or the warranty will be invalid.

- ▲ **Only use batteries sizes #13, #313, #675 Ni-Mh or equivalent size Ni-Mh rechargeable batteries.**

- ▲ **Do not disassemble.**

This Battery charger contains no serviceable parts.

- ▲ **Observe cleaning instructions.**

Do not clean with solvents such as Benzine, Thinners, Turpentine, Methylate Spirits etc. Use a soft dry cloth to wipe clean.

- ▲ **Do not immerse in water**

This product contains electrical components which can be damaged if they are exposed to water or high levels of humidity.

▲ Do not incinerate.

This product contains electrical components which can explode if exposed to flame or intense heat.

▲ Place in direct sunlight to recharge.

It is suggested that hearing aids and rechargeable hearing aid batteries be charged at night from internal batteries and that the internal AA rechargeable batteries be charged in the day.

▲ If charging using household light, insure the solar charger is within 30cm or 12 inches from the light source.

▲ Only government approved (examples ANSI, CE, CSA) electrical wires must be used if charging electrically. Most but not all approved cell phone connections will work.

▲ During charging the LED will illuminate.

When inserting the hearing aid and/or rechargeable batteries the green LED will be on during recharge. When the green LED is off it means that battery charge is complete. If the red LED is on it means its polarity was inverted.

▲ Keep out of reach of small children, if batteries are swallowed seek medical attention immediately.

Using your Solar Charger

Inserting Batteries

Step 1 - Insert batteries with (+,-) as shown.
(See fig. 1)


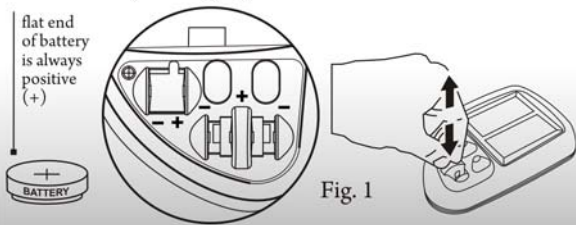
Fig. 1


Inserting Solar Ear Rechargeable Hearing Aid

Step 2 - Insert hearing aid as shown, make sure the (+,-) is the same as in Fig1.

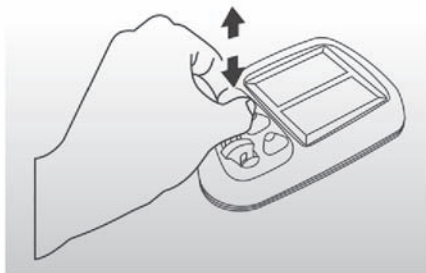Removing Batteries

Step 3 - Push button down to eject batteries from recess, the refer to (Step 1) for removal.

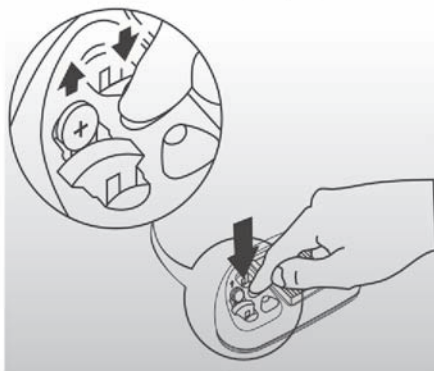Warranty

Solar Ear Charger

Limited Warranty - 2 years

The manufacturer warrants that this product is free, under normal and maintenance, from defects in material and workmanship. If any such defects should be found in this product within the warranty period the manufacturer shall, at its option, repair or replace the product as specified below:

The warranty shall not apply to:

- a. Any defects caused or repairs required as a result of abusive operation, negligence, accident, mishandling or inappropriate use as outlined in the operation instructions.
- b. Any Solar Ear tampered with, modified, adjusted or repaired by any party other than the manufacturer or authorized distributor.
- c. Damage caused or repairs required as a result of the use with items not specified or approved by the manufacturers.
- d. Any product received without appropriate serial number.
- e. Defects caused by battery leakage.

**THIS CARD MUST ACCOMPANY YOUR
INSTRUMENT
FOR WARRANTY SERVICE**

Should Solar Charger fail to operate during the warranty period, warranty service may be obtained upon delivery of Solar Ear direct to the manufacturer or authorized distributor. The sender will be responsible for any transportation and insurance costs incurred.

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Manufactured and assembled by young adults who are deaf.

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