

## TECHNICAL DATA

# BW-265




## Description

BW-265 is the most attractive choice in Body Worn hearing aids. Adjustable volume controls for each ear deliver clear and balanced sound and great flexibility for special needs. The audio input function makes it possible to connect the B80 with external audio devices. Enjoy the simple life of being reconnected to the world of sound.

## What's Special?

### BW-265

- \* Two channels
- \* Individual-ear adjustment
- \* Binaural headset
- \* New Class- AB Amplifier
- \* Low Harmonic Distort
- \* Low noise
- \* Audio input support
- \* Telecoil
- \* With trimmers
- \* Powerful Output
- \* MTO Switch

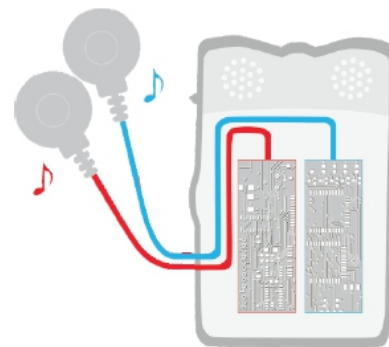


The dual amplification design enables users to adjust the volume to each ear independently for both balance and special circumstances. Now you just need one hearing instrument for your two ears.

## TECHNICAL DATA


# BW-265



### Fitting range

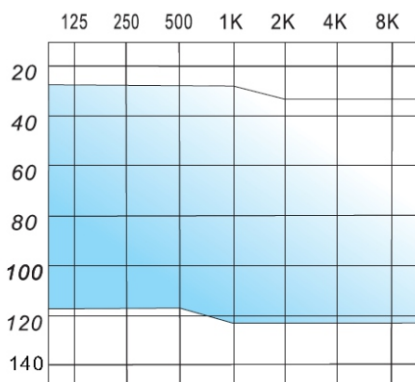

### Testing diagram

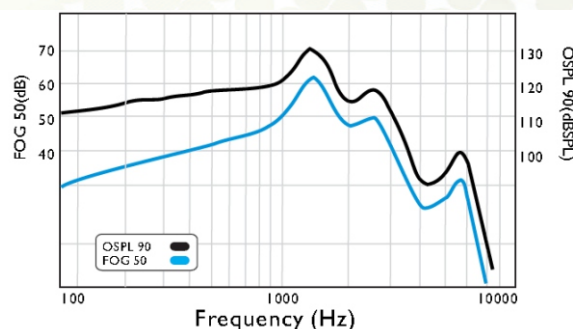

### Main Technical Data

Technical Data	BW-265
IEC-118-7/2005(2cc-coupler)	
Peak OSPL 90(db SPL)	132
HAF OSPL 90(db SPL)	122
Peak Gain(db)	65
HAF full on Gain(db)	55
Total Harmonic Distort:	
500Hz(%)	< 3.8
800Hz(%)	< 2.5
1600Hz(%)	< 1.0
EQ Input Noise(db)	< 25
Battery Current(mA)	< 4.5
Signal Direction	Analog
T-Coil	Y
Potentiometer Features	VC
Trimmers	H+P+H+P+AU
Cord	Binaural headset
Number of Channels	
Battery type	AAA
Module dimensions	660.x43.1x20.3mm
Colors Available	Black-Silver

The measurements of the represented curves are in default setting of the device. All trimmers and potentiometers are set in full-on position. All technical(typical) data are subject to change without prior notice.